## **Hydrazide-Terminated Magnetic Beads**

#### Introduction

BcMag™ Hydrazide-Terminated Magnetic Beads are uniform magnetic beads grafted with a high density of hydrazide functional groups on the surface. Hydrazide chemistry is effective for labeling, immobilizing, or conjugating glycoproteins via glycosylation sites, which are frequently (as with most polyclonal antibodies) positioned in domains distant from the critical binding sites whose function needs to be preserved. Coupling antibodies in this manner selectively targets heavy chains in the Fc portion of the molecule, assisting in the best possible preservation of antigen binding activity by the ends of the Fv regions. At pH 5 to 7, hydrazide-terminated supports and compounds will conjugate with carbonyls of oxidized carbohydrates (sugars), resulting in the production of hydrazone linkages (Fig.1). The hydrazide beads are suitable for conjugating larger glycoproteins and glycolipids, carbohydrates, or other ligands. Moreover, the hydrophilic surface ensures low nonspecific adsorption, excellent dispersion, and easy handling in various buffers.

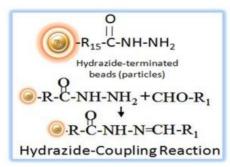


Fig.1 Hydrazide coupling reaction.

#### Workflow

The Beads work perfectly as affinity resin for affinity purification to refine molecules, cells, and parts of cells into purified fractions. After conjugation with ligands, add the beads to a sample containing the target molecules, then mix, incubate, wash and elute the target molecules (Fig.2)



#### Features and benefits.

- Specific immobilization—the hydrazide-activated beads bind exclusively to pure glycoproteins containing sugar groups that have been gently oxidized with periodate (e.g., sialic acid).
- Stable covalent bond with low levels of ligand leakage
- Maintains antibody function—immobilizes IgG via the Fc region, leaving both antigen binding sites available for target capture.
- Low nonspecific binding
- High capacity—Immobilize 15-20µg antibody/mg beads.

Specification		
Composition	Magnetic Bead grafted with hydrazide group on the surface.	
	Short Term (<1 hour): pH 3-11; Long-Term: pH 4-10	
Stability	Temperature: 4°C -140°C; Most organic solvents	
Magnetization	~40-45 EMU/g	
Type of Magnetization	Superparamagnetic	
Formulation	Lyophilized Powder	



#### **Instruction Manual**

Functional Group Density	1µm Magnetic Beads	~200 µmole / g of Beads
	5μm <b>Magnetic</b> Beads	~180 µmole / g of Beads
Storage	Ship at room temperature. Store at 4°C upon receipt. Do not freeze.	

#### Protocol Number of Beads

#### Note:

- This protocol can be scaled up as needed. We strongly recommended titration to optimize the number of beads used for each
  application.
- Avoid tris or other buffers containing primary amines because these will compete with the intended coupling reaction.

#### **Materials Required**

- 1. Magnetic rack (for manual operation): Based on sample volume, the user can choose one of the following magnetic Racks: BcMag rack-2 for holding two individual 1.5 ml centrifuge tubes (Cat. # MS-01); BcMag rack-6 for holding six individual 1.5 ml centrifuge tubes (Cat. # MS-02); BcMag rack-24 for holding twenty-four individual 1.5-2.0 ml centrifuge tubes (Cat. # MS-03); BcMag rack-50 for holding one 50 ml centrifuge tube, one 15 ml centrifuge tube, and four individual 1.5 ml centrifuge tubes (Cat. # MS-04); BcMag<sup>TM</sup> rack-96 for holding a 96 ELISA plate or PCR plate (Cat. # MS-05).
- 2. Coupling Buffer: 0.1 M sodium phosphate, pH 7.0
- 3. Oxidizing Agent: Sodium meta-Periodate (NaIO4) (Sigma, Cat# S1878)
- 4. Washing Buffer: 1M NaCl
- 5. PBS

#### **Ligand Coupling**

#### A. Beads preparation

- 1. Weight 30 mg beads and add them into a 1.5 ml centrifuge tube.
- 2. Add 1 ml coupling buffer and resuspend the beads very well by vortexing or pipetting.
- 3. Insert the tube into a magnetic rack for 1-3 minutes until the supernatant becomes clear. Aspirate and discard the supernatant with a pipette while the tube remains in the rack.
- 4. The beads are ready for coupling.

#### B. Oxidation of Glycoprotein or other ligands

**Note:** The reaction is light sensitive and should be performed in the dark.

- 1. Dissolve or dilute 0.5-10 mg glycoprotein or other ligands in 1 ml coupling buffer. (Note: If the protein or other ligands is already suspended in other buffers, perform a buffer exchange by dialysis or a desalting column.)
- Add the protein or other ligands solution to an amber vial containing 2 mg sodium meta-periodate (final concentration10mM). Swirl gently to dissolve the oxidizing agent.)
- 3. Incubate the sample in the dark at room temperature for 45 minutes with good mixing (end-over-end).

#### C. Conjugation

 Add the oxidized protein or other ligands solution to the prepared magnetic beads from Step A4 and mix well by vortexing or pipetting.

**Note:** Coupling efficiency depends on the structure and the size of the target glycoprotein or other ligands. The user should empirically optimize the ratio of the protein to the beads.

- 2. Incubate the reaction **in the dark** at room temperature overnight with good mixing (en-over-end).
- 3. Insert the tube into a magnetic rack for 1-3 minutes until the supernatant becomes clear. Aspirate and discard the supernatant with a pipette while the tube remains in the rack. Remove the tube from the rack and resuspend the beads with 5 ml Coupling Buffer by vortex or pipette.
- 4. Repeat steps 3 for three times.
- 5. Resuspend the beads in PBS buffer with 0.01% azide (w/v) to desired concentration and store at 4°C until use. Do not freeze.



#### **Instruction Manual**

#### D. General affinity purification Protocol

#### Note:

- This protocol is a general affinity purification procedure. Designing a universal protocol for all protein purification is impossible because no two proteins are precisely alike. To obtain the best results, the user should determine the optimal working conditions for purifying the individual target protein.
- We strongly recommended titration to optimize the number of beads used for each application based on the amount of the target protein in the crude sample. Too many magnetic beads used will cause higher backgrounds, while too few beads used will cause lower yields. Each mg of magnetic beads typically binds to 10-20 µg of the target protein.
- 1. Transfer the optimal amount of the beads to a centrifuge tube. Place the tube on the magnetic rack for 1-3 minutes. Remove the supernatant while the tube remains on the rack.
- 2. Remove the tube and wash the beads with 5-bed volumes of PBS buffer by vortex for 30 seconds. Leave the tube at room temperature for 1-3 minutes. Place the tube on the magnetic rack for 1-3 minutes. Remove the supernatant while the tube remains on the rack.
- 3. Repeat step 2 two times.
- 4. Add washed beads to the crude sample containing the target protein and incubate at room or desired temperature for 1-2 hours (Lower temperatures require longer incubation time).
  - **Note:** Strongly recommended to perform a titration to optimize incubation time. More prolonged incubation may cause higher background.
- 5. Extensively wash the beads with 5-beads volumes of PBS buffer or 1M NaCl until the absorbance of eluting at 280 nm approaches the background level (OD 280 < 0.05).
  - Note: Adding a higher concentration of salts, nonionic detergent, and reducing agents may reduce the nonspecific background. For example, adding NaCl (up to 1-1.5 M), 0.1-0.5% nonionic detergents such as Triton X 100 or Tween 20, and a reducing reagent such as DTT or TCEP (we usually use 3mM) to the washing buffer.
- 6. Elute the target protein by appropriate methods such as low pH (2-4), high pH (10-12), high salt, high temperature, affinity elution, or boiling in an SDS-PAGE sample buffer.



# **Instruction Manual**

Related Products		
Amine-Terminated Magnetic Beads	Iodoacetyl-Activated Magnetic Beads	
DADPA-Activated Magnetic Beads	Peptide conjugation buffer Kit-I	
Carboxyl-Terminated Magnetic Beads	Peptide conjugation buffer Kit-II	
Epoxy-Activated Magnetic Beads	DVS-Activated Magnetic Beads	
Hydrazide-Terminated Magnetic Beads	NHS-Activated Magnetic Beads	
Glycoprotein and Antibody Conjugation Kit-I	Hydroxyl-Terminated Magnetic Beads	
Glycoprotein and Antibody Conjugation Kit-II	Sulfhydryl-Terminated Magnetic Beads	
Aldehyde-Activated Magnetic Beads	Tosyl-Activated Magnetic Beads	
Silica-Modified Magnetic Beads	CDI-Activated Magnetic Beads	
Alkyne-Activated Magnetic Beads	Thiol-Activated Magnetic Beads	
Azide-Activated Magnetic Beads	Cleavable NHS-Activated Magnetic Beads	
Cleavable Amine-Terminated Magnetic Beads	Cleavable Azide-Activated Magnetic Beads	
Cleavable Carboxyl-Terminated Magnetic Beads	Cleavable Alkyne-Activated Magnetic Beads	
Cleavable Epoxy-Activated Magnetic Beads	Cleavable Iodoacetyl-Activated Magnetic Beads	
Cleavable Hydrazide-Terminated Magnetic Beads	Cleavable Tosyl-Activated-Magnetic Beads	
Cleavable Aldehyde-Activated Magnetic Beads	Streptavidin Magnetic Beads	
Boronate Affinity Magnetic Beads	Cleavable Streptavidin Magnetic Beads	
Monomer Avidin Magnetic Beads		